

USB® VeriQuest™ Fast Probe qPCR Master Mix (2X)

Applications:

- Amplification from genomic DNA or cDNA input
- Gene expression validation

Features:

- Uses fast mode thermal cycling conditions for results in one-third of the time when compared to standard run protocols
- One-tube master mix with ROX™ Reference Dye included—just add template, primers and water
- Reproducibility and consistency over a broad dynamic range: 8 orders of magnitude linear detection range (Fig. 1)
- High specificity with exceptional performance on challenging GC-rich regions
- High sensitivity and precision with limited targets (Fig. 2)
- Stable at room temperature for 72 hours in a pre-assembled reaction (Fig. 3)
- Contains dUTP and Uracil-DNA Glycosylase (UDG or UNG) for carry-over contamination prevention
- Compatible with a wide variety of real-time PCR instruments


USB VeriQuest Fast Probe qPCR Master Mix is optimized for SYBR® Green detection on all instruments that utilize ROX as a passive reference dye. This master mix can be run using Fast mode cycling protocols, providing results in one-third of the time when compared to standard real-time PCR reagents. The 2X ready-to-use master mix contains hot start VeriQuest Taq DNA Polymerase, MgCl₂, ultrapure nucleotides with an optimized dUTP:dTTP ratio, Uracil-DNA Glycosylase (UDG or UNG), and ROX Passive Reference Dye in a proprietary reaction buffer. The hot start Taq polymerase has no polymerase activity prior to the initial heat activation step which allows reaction assembly at room temperature. The mix offers higher specificity and sensitivity, virtually eliminating non-specific primer amplification which can negatively affect the efficiency and accuracy of the data (Fig. 1). Since the mix contains dUTP and UDG, carry-over contamination prevention can be performed prior to amplification. VeriQuest Fast Probe qPCR Master Mix offers the same sensitivity and specificity found with our standard mode mixes with results in a fraction of the time.

VeriQuest Fast Probe qPCR Master Mix shows exceptional efficiency and high specificity on challenging templates such as high GC- and AT-rich regions. The master mix is also capable of duplex target detection without compromising the integrity of the data (Fig. 4). The optimized formulation ensures no sacrifice in quality with increased speed.

The sensitivity of the master mix allows for discrimination from a 1.33-fold difference in gene target amount detected. In Figure 2, a 1.33 to 10-fold dilution series of 10 ng to 1 ng of cDNA reverse-transcribed from HeLa total RNA were amplified with an efficiency of >98.8%.

VeriQuest Fast Probe qPCR Master Mix is highly stable and easy to work with. It remains stable at room temperature for 72 hours in a pre-assembled reaction and can be stored at 4°C for convenient handling. The mix also allows for room temperature reaction set-up. The speed and stability make VeriQuest Fast Probe qPCR Master Mix ideal for high-throughput handling.

Fig. 1. Linear detection range

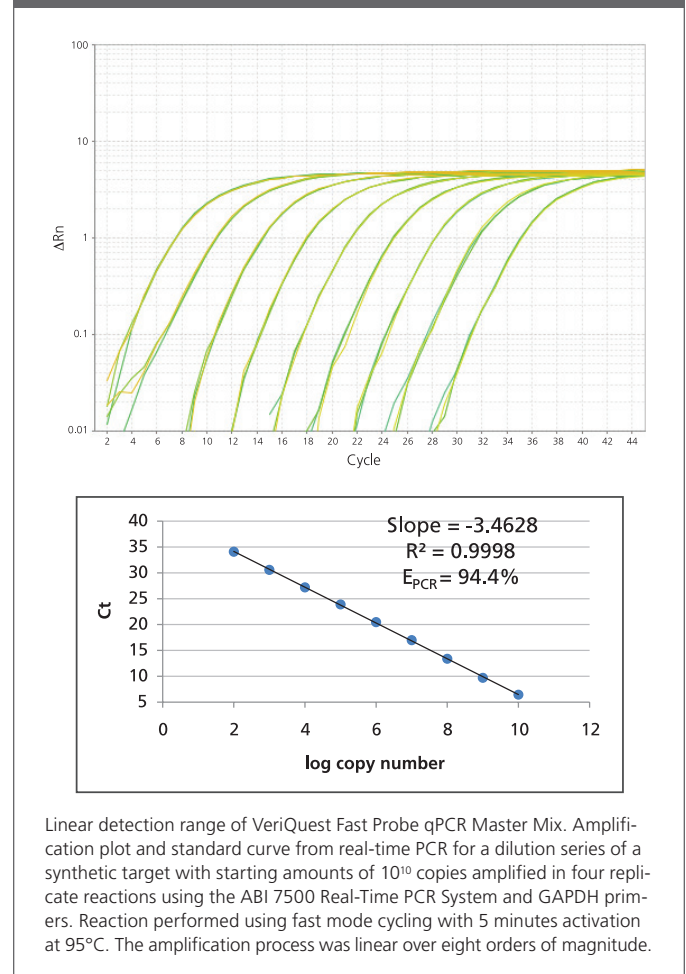

Linear detection range of VeriQuest Fast Probe qPCR Master Mix. Amplification plot and standard curve from real-time PCR for a dilution series of a synthetic target with starting amounts of 10¹⁰ copies amplified in four replicate reactions using the ABI 7500 Real-Time PCR System and GAPDH primers. Reaction performed using fast mode cycling with 5 minutes activation at 95°C. The amplification process was linear over eight orders of magnitude.


Components:

75680	100 reactions	1 ml
	500 reactions	5 ml
	1,000 reactions	2 x 5 ml
	2,500 reactions	5 x 5 ml
	5,000 reactions	10 x 5 ml

VeriQuest Fast Probe qPCR Master Mix (2X)


Product Code	Pack Size
75680	100 reactions
	500 reactions
	1,000 reactions
	2,500 reactions
	5,000 reactions

Fig. 3. VeriQuest Fast Probe qPCR Master Mix stability

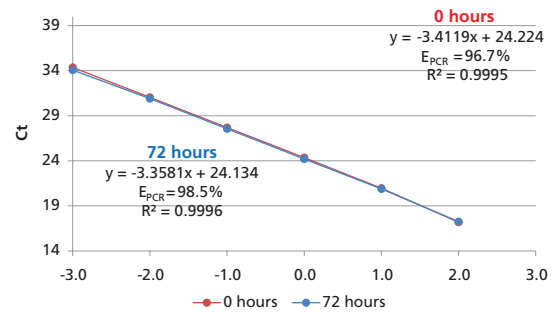


GAPDH was detected from pre-assembled PCR reactions incubated at room temperature for 72 hours. Results were compared to the freshly prepared mix (0 hours).

Fig. 2. High sensitivity and precision

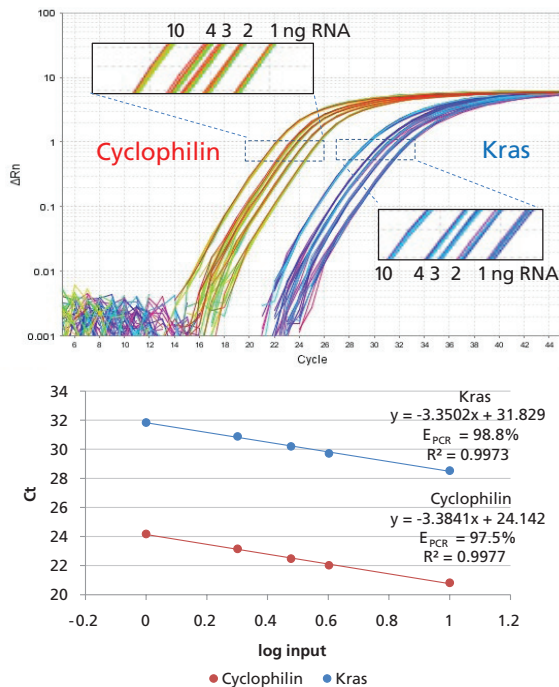

High sensitivity and precision in limited target quantification. Amplification plot (bottom) and standard curve (top) from real-time PCR for a 1.33 to 10-fold dilution series of 10 ng to 1 ng of cDNA reverse-transcribed from HeLa total RNA.

Fig. 4. Single-target and duplex PCR amplification

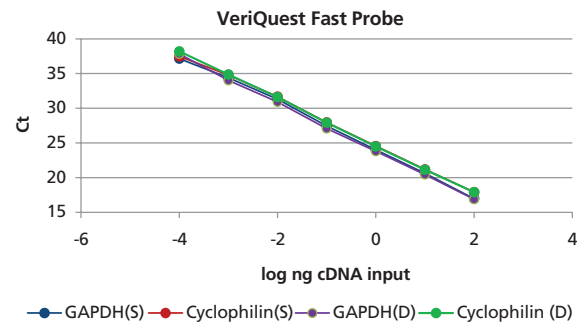

USB VeriQuest Fast Probe Mix

	slope	E _{PCR}	R ²
GAPDH(S)	-3.4039	96.7%	0.9983
Cyclophilin(S)	-3.3229	100.0%	0.9986
GAPDH(D)	-3.4690	94.2%	0.9996
Cyclophilin(D)	-3.4069	96.6%	0.9998

Fast cycling conditions

USB VeriQuest Probe Mix

	slope	E _{PCR}	R ²
GAPDH(S)	-3.4507	94.9%	0.9996
Cyclophilin(S)	-3.4116	96.4%	0.9999
GAPDH(D)	-3.3623	98.3%	1.0000
Cyclophilin(D)	-3.3996	96.9%	1.0000

Standard cycling conditions

VeriQuest Fast Probe qPCR Master Mix single-target and duplex PCR amplification of Cyclophilin Fam-BHQ-1® probe and GAPDH Cy5-BHQ-1 probe. Single-target (S) and duplex reactions (D).

For research use only. Not for use in diagnostic procedures.

Affymetrix and USB are registered trademarks of Affymetrix Inc. VeriQuest is a trademark of Affymetrix Inc. Taq DNA Polymerase—sold under licensing arrangements with Applied Biosystems. Purchase is accompanied by a limited license to use it in the Polymerase Chain Reaction (PCR) process in conjunction with a thermal cycler whose use in the automated performance of the PCR process is covered by the up-front license fee, either by payment to Perkin-Elmer or as purchased, i.e., an authorized thermal cycler. ROX is a trademark of Applied Biosystems or its subsidiaries in the US and certain other countries. SYBR is a registered trademark of Molecular Probes, Inc. and is provided under an agreement with Molecular Probes, Inc. BHQ-1 is a registered trademark of Biosearch Technologies.

Affymetrix, Inc.

USB® Products
 26111 Miles Road
 Cleveland, Ohio 44128
 Tel: 800-321-9322 | 216-765-5000
 Fax: 800-535-0898 | 216-464-5075
 USBcustomerserv@affymetrix.com
usb.affymetrix.com

Affymetrix UK Ltd.

USB® Products
 Voyager, Mercury Park,
 Wycombe Lane, Wooburn Green,
 High Wycombe HP10 0HH, United Kingdom
 Tel: +44 (0)1628 55 2600
 Fax: +44 (0)1628 55 2675
 USBcustomerserveurope@affymetrix.com

Affymetrix Pte Ltd

USB® Products
 7 Gul Circle, #2M-01
 Keppel Logistics Building
 Singapore 629563
 Tel: +65 63957301
 Fax: +65 63957300
 USBcustomerserv@affymetrix.com